

Vivera Monobloc M301

MEDIROL

The powered monobloc represents a compact and simple solution for the daily transport of ordinary and bariatric patients.

Key features and benefits

- Battery-powered hydraulic lift system
- Integrated battery with battery management system
- Automatic battery charging in fastener system
- Settable load height
- High visibility powder-coated frame
- Integrated bumper system
- Retractable head section
- One-hand release, infinite positioning, pneumatically assisted backrest and legrest
- Rubber wheels for better shock absorption
- Steer lock system enables wheels either to swivel or to lock into straight-line position
- Foot-end wheel locks



320 kg

MAXIMUM LOAD CAPACITY



5 Ah

BATTERY CAPACITY



Up to 95 cm

LOADING HEIGHT




Specifications

Type	M301	Articulation	
Length		Backrest	0° - 90°
Standard	200 cm	Shock position	0° - 30°
Minimum	160 cm	Headrest ²	± 35°
Width	60 cm	Kneerest ²	0° - 50°
Height		Head extension ²	+29 cm
Highest position	133 cm	Extended Patient Surface ²	94 cm
Lowest position	38 cm	Battery capacity	5 Ah
Weight ¹	75 kg	Maximum load capacity	320 kg
Wheels		Recommended loading height	< 95 cm
Diameter	15 cm	Recommended fastener system ³	F401
Width	5 cm		

¹ Stretcher is weighed with battery pack, without mattress and restraints. ² Optional. ³ Meets dynamic crash standards EN-1789 with F401 fastening system. Specifications are rounded to the nearest whole number. MEDIROL reserves the right to change specifications without notice.

IP 66 EN 1789 EN 1865 CE

Vivera Fastener F401

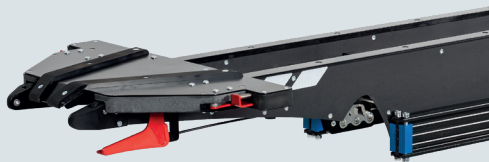

VISIT [MEDIROL.CZ](https://www.medirol.cz) TO VIEW THE ENTIRE MEDIROL ACCESSORIES

